


FIBER FORCE®

Light cure glass fiber & resin
for the dental laboratory



Metal-free fiber reinforcement for:

- fabricating new dentures and repairing dentures
- reinforcing acrylic elements
- reducing denture fracture risk

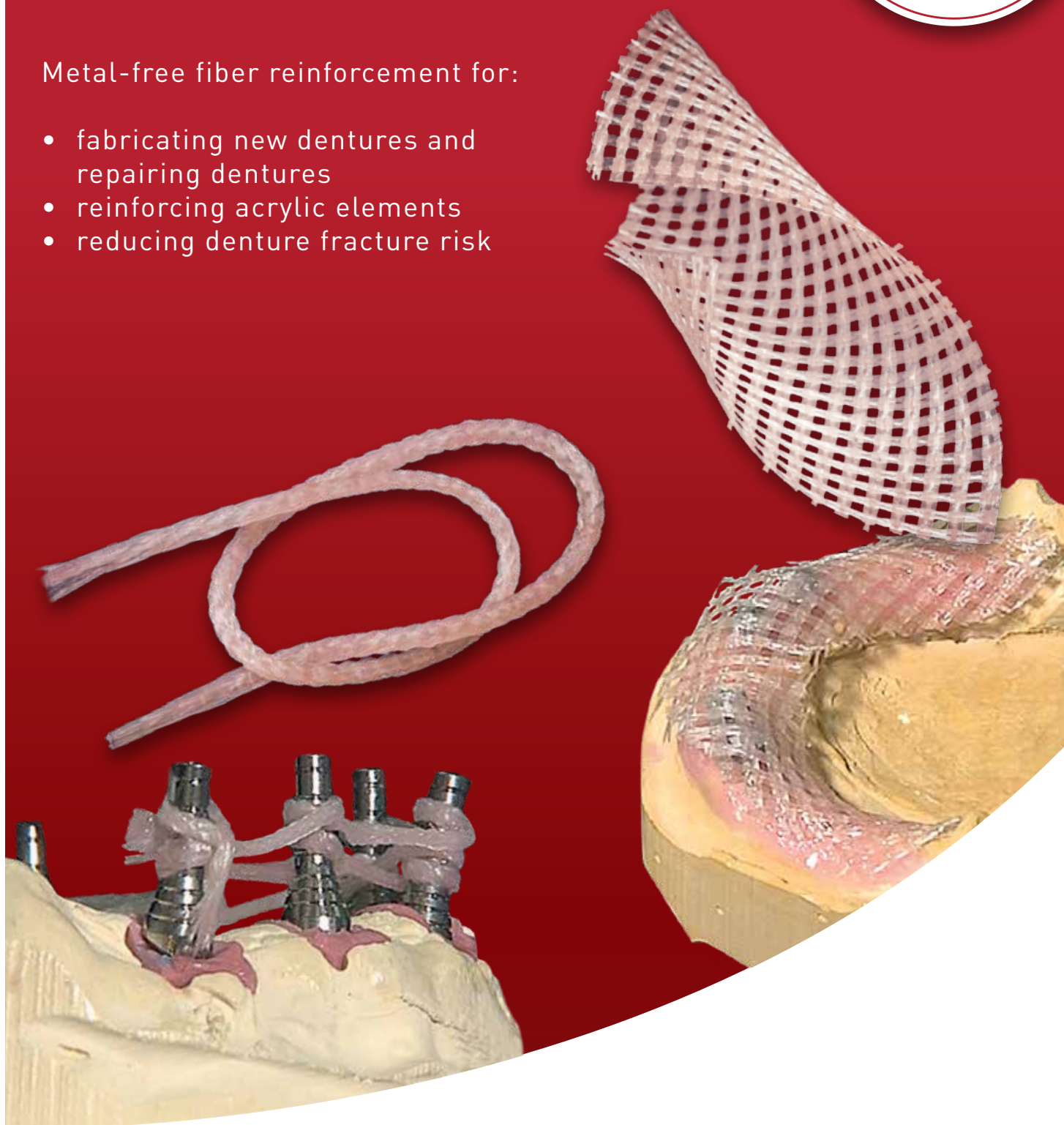

FIBER FORCE®


Fabricate stable, metal-free frameworks for screw-retained implant-supported dentures



**Pink, metal-free
Light Cure Braids
(2 pc.) Ø 2.2 mm x 150 mm**

REF CI08002

- For denture repairs
- Supports implant-supported dentures
- Reinforces acrylic dentures

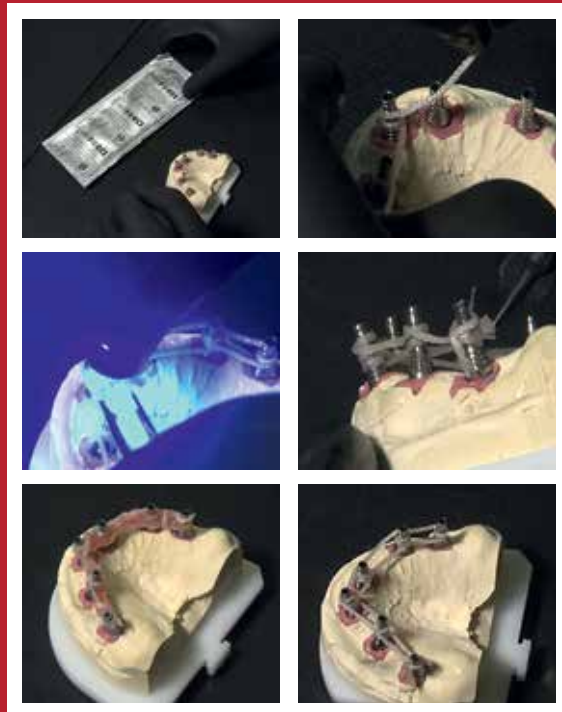


Reinforce removable dentures



**Pink, metal-free
Light Cure Mesh
(2 pc.) 80 mm x 53 mm**

REF CI08004

- For new full dentures
- For implant-supported dentures
- For denture repairs

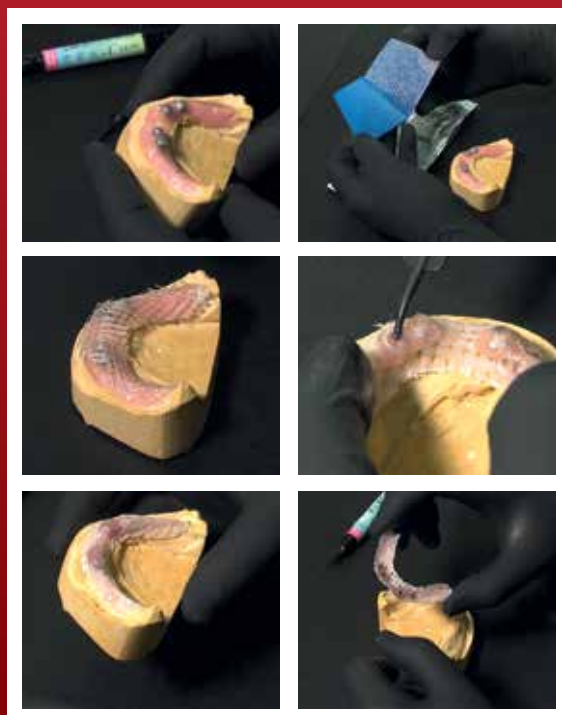



€ 166*

EZ VAC Vacuum Unit

REF CI13002

The EZ VAC is an easy-to-use vacuum forming device specifically designed for FIBER FORCE® mesh. EZ VAC allows you to fabricate perfectly fitted, thin FIBER FORCE® reinforcements anytime. A silicone membrane pulls down the uncured mesh, which can then be finely adjusted to the model. Then, the EZ VAC is inserted into a light-curing unit, and the mesh is cured into its final shape.



€ 326*

Kit Overdenture/Partial Denture

REF CI17001

- 1 x EZ VAC Vacuum Unit
- 2 x Pink Light Cure Braid
- 2 x Pink Light Cure Mesh
- 1 x Pink Light Cure Resin



Light Cure Resin 3 ml-syringe

- For fabricating mesh
- For tacking the framework to the model
- For repairing as a replacement for acrylic

pink
REF CI08027

tooth-colored
REF CI08028

€ 28*
each



€ 221*

CST Starter Kit

REF CI08111

Self-supporting three-dimensional fiber structure for implant-supported bridge

- 1 x Pink Hybrid Fiber 1:6 (Ø 1.3 mm x 450 mm)
- 1 x Pink Hybrid Fiber 1:4 (Ø 1.2 mm x 300 mm)
- 3 x Support Pillars
- 1 x Bond "James-2"
- 1 x Pink Light Cure Resin (syringe)

FIBER FORCE® is a registered trademark of Bio Composants Médicaux.

* Prices + VAT + shipping

You can find new products here:

[scan QR-Code]



Order online in the Champions-Implants-Shop or per Fax!

champions  implants

Order

FIBER FORCE®

_____ piece(s) Pink Light Cure Braids (REF CI08002)	€ 43/pc.*
_____ piece(s) Pink Light Cure Mesh (REF CI08004)	€ 59/pc.*
_____ piece(s) Pink Light Cure Resin (REF CI08027)	€ 28/pc.*
_____ piece(s) White Light Cure Resin (REF CI08028)	€ 28/pc.*
_____ piece(s) CST Starter Kit (REF CI08111)	€ 221/pc.*
_____ piece(s) EZ VAC Vacuum Unit (REF CI13002)	€ 166/pc.*
_____ piece(s) Kit Overdenture/Partial Denture (REF CI17001)	€ 326/pc.*

Billing address (eventually with stamp)

Name _____
Street _____
ZIP Code/City/Country _____
Phone _____ Fax _____
Customer No. _____

Shipping address (eventually with stamp)

Name _____
Street _____
ZIP Code/City/Country _____
Phone _____ Fax _____
Customer No. _____

* Prices + VAT + shipping

Fax-Order: 0049 6734 10 53

Free shipping, EU-Countries only, for orders from € 500 (net price) up.

By placing an order for the following references, the customer accepts the General Terms and Conditions (GTC) of Champions-Implants GmbH (GTC on the homepage > Shop > GTC).

Please send me the current product catalog

CHAMPIONS-IMPLANTS GMBH

Im Baumfeld 30 | Champions Platz 1 | D-55237 Flonheim
fon +49 (0)6734 91 40 80 | fax +49 (0)6734 10 53
info@champions-implants.com

Join the courses offered by the Future Dental Academy!

- MIMI®-MotivaCtion
- Prosthodontics & Orthodontic Aligners
- Biological Implantology

You can find dates and registration forms on our website:

www.champions-implants.com



Bestellschein 4-Seiter Fiber Force 2020 EN
16. April 2020, 11:34 vorm.

BE PART OF THE CHAMPIONS® (R)EVOLUTION!

More information about our products:

www.champions-implants.com

champions  implants